

## Yellofins® Elite Stirrups


Supports up to 500 lbs (227 kg)




Easily Grip and Squeeze...



...Raise and Abduct...

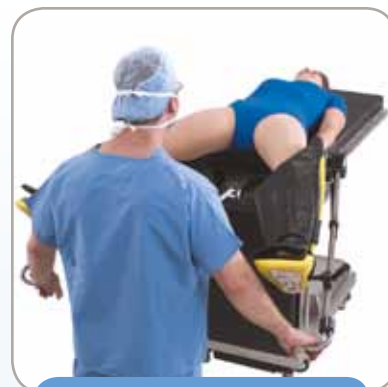

...to Any Position.

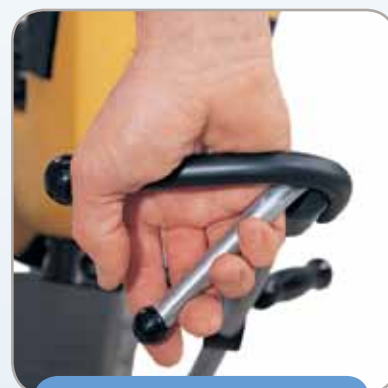

Squeeze and Position

## Yellofins® Stirrups

Supports up to 350 lbs (159 kg)

## Simultaneous Lithotomy and Abduction

**Note: Stirrups are Sold in Pairs**








The Yellofins Stirrup represents the gold standard in lithotomy positioning. Yellofins Stirrups allow for easy adjustment of abduction and lithotomy while maintaining the sterile field. The boot is extended along its lateral side (referred to as a 'lateral fin') and is padded on the interior. The purpose of this fin is to cover and protect the head of the fibula and hence, the peroneal nerve. The Yellofins Stirrup enables safe and easy positioning while providing enhanced surgical site access.

- **Boot Design:** Extended lateral fin is designed to cover and protect the head of the fibula and hence, the peroneal nerve.
- **Floating Boot:** Self-adjusting boot minimizes pressure on the calf, when moving the stirrup.
- **Lift-Assist™ Technology:** Allows easy movement of the leg when placing it in the desired position.
- **Yellofins Reusable Pad:** Completely encapsulates the foot, ankle and calf for comfortable and secure positioning.
- **Squeeze Grip Handle:** Provides easy intraoperative adjustment without compromising the sterile field. Simply release the handle to secure the leg holder in all directions.
- **Lithotomy Range:** Yellofins, and Yellofins Elite with Lift-Assist technology – Set the stirrup in any position between +84° to –33° lithotomy. Yellofins without Lift-Assist technology – Set the stirrup in any position between +90° to –35° lithotomy.
- **Abduction Range:** Set the stirrup in any position between +25° to –9° abduction.
- **Lithotomy & Length Indicators:** Visual indicators that help ensure precise positioning.

#O-YFES Yellofins Elite Stirrups with Lift-Assist (500 lbs/227 kg Capacity)

#O-YFASI Yellofins Stirrups with Lift-Assist (350 lbs/159 kg Capacity)

#O-YFSI Yellofins Stirrups without Lift-Assist (350 lbs/159 kg Capacity)

What's Included	 <p>Yellofins Elite Stirrups with Lift-Assist #O-YFES 1 Pair</p>	OR	 <p>Yellofins Stirrups with Lift-Assist #O-YFASI 1 Pair</p>	OR	 <p>Yellofins Stirrups without Lift-Assist #O-YFSI 1 Pair</p>
What's Required	 <p>Yellofins Reusable Pad #O-AMA68S0.6S 1 Pair Required</p>		 <p>Easy Lock Blade Clamp #A-40040 (US), #A-40041 (EU), #A-40042 (UK), #A-40043 (Jpn) #A-40044 (Den), 2 Required Specify Rail Type</p>		 <p>Safety Drape #A-22000 Case of 10 Pairs</p>
What's Optional	 <p>Easy Catcher Disposable System (Case of 10) (Sterile) (Non-Sterile) (See Page 85)</p>		 <p>Stirrup Cart #A-30015</p>		 <p>Extended Warranty #A-XWARYFPRO Available in US Only</p>